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massage to go

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BellaBambi consists of biocompatible silicone, making it dermatologically safe and suitable also for allergy sufferers. BellaBambi is made from natural quartz sand and can be recycled. //

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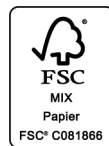
Thank you for buying BellaBambi

Congratulations on buying BellaBambi. This is a multi-purpose and durable product for beauty, relaxation and sports. BellaBambi is dermatologically safe, consists of biocompatible silicone and is produced exclusively in Germany. To have as much success and lasting pleasure with your product, we kindly ask you to read this instructions carefully.

Thank you!
Your team from BellaBambi

BellaBambi increases blood circulation, improves lymphatic circulation and promotes flexibility. It raises the different layers of the skin, stretches fascia, can release adhesions and stimulates the metabolism. As you drag BellaBambi over your skin you can also achieve an exfoliating effect.

Recommendations

For most applications, an oil to assist lubricity is helpful. BellaBambi can also be used with hyaluronfluid, cellulite creams or body lotions. Many suitable products can be found on www.bellabambi.com

Please note our experiences with possible skin reactions: when BellaBambi is permanently left in place:

	pressure mark	redness	effusion
sensitive skin	approx. 1-2 min	approx. 2-3 min	approx. 5 min
normal skin	approx. 2 min	approx. 4 min	approx. 7 min
robust skin	approx. 3 min	approx. 6 min	approx. 9 min

Bamberger Wellness GmbH
Industriestr. 25
D-89423 Gundelfingen
Germany
+49 9073 9210461
info@bellabambi.de
www.bellabambi.de

Precautions / Contraindications

Do not apply BellaBambi on nipples, sexual organs, anus, eyes, eyelids, lips, mouth, in the axilla and in the inguinal region.

Please consult with your physician before using BellaBambi if one or more of the following conditions apply: Pregnancy, acute injuries, open wounds, inflammations, varicose veins, disposition to hematomas and ulcers caused by diabetes or vascular diseases, irritable skin disorders such as psoriasis or urticaria (hives), with blood-thinning medication therapies, tumors, Complex regional pain syndrome (CRPS), as well as any unclear or acute medical conditions.

Intensities

SENSITIVE for soft application
VITALITY for vital application
INTENSE for strong application

Cleaning and Disinfection

BellaBambi can be cleaned very easily in the dishwasher or with mild soap suds. After drying, please treat BellaBambi with a suitable disinfectant for silicones (without phenol derivatives or alkyl amines (Glucoprotamin). We recommend: Desomed® Rapid AF). Suitable cleaning agents and disinfectants can be found on www.bellabambi.com

Internet / Social Media

www.bellabambi.de
facebook.com/bellabambi.de
pinterest.com/BambergWellness
youtube.com/bambergervacu
beautyrelaxandsports.tumblr.com

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Muscle oil with skin caring properties. Relaxes cramped muscles and supports the BellaBambi treatment.

Ingredients: Hypericum perforatum oil, Glycine soja oil, Arnica montana flower, Rosmarinus officinalis oil, Juniperus communis oil, Eucalyptus globulus oil, Pinus pumilio oil, Lavendula angustifolia oil, Citral*, Citronellol*, Coumarin*, Farnesol*, Limonene*, Linalool*, Cinnamal*, Isoeugenol*.) from natural essential oils

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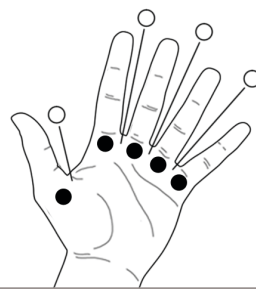
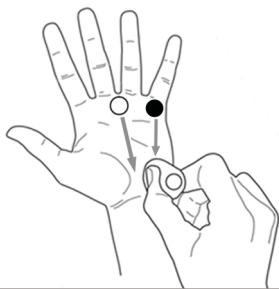
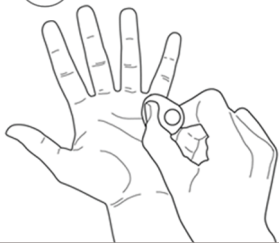
Place BellaBambi on given marks and drag it to your wrist.

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Place BellaBambi in the middle of your palm. Flex your fingers, almost to a fist, and extend them again. Please do this very slowly.

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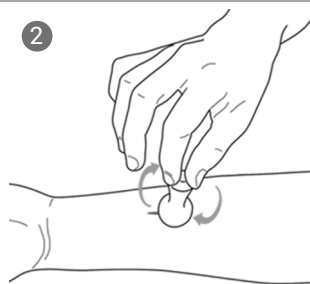
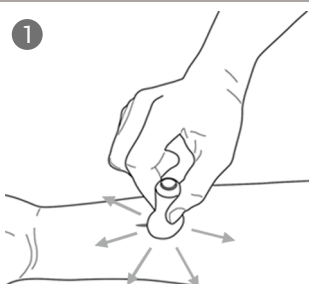
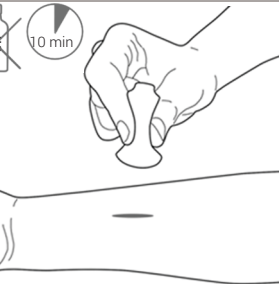
1 Move BB over your scar in a "star" pattern.
2 Place BB on top of your scar and describe circles with the tip of BellaBambi, or push it into one direction for about 20 seconds.

Important: Please consult your doctor first and do not start before your scar is clean, dry and closed.

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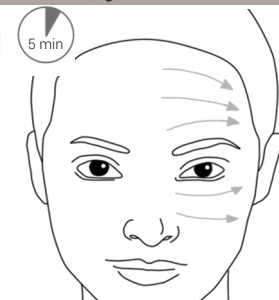
Place BellaBambi below your hairline in the middle of your forehead and drag it outward. Relocate BellaBambi 1 cm below and repeat until you reach your eyebrow. Move to your cheek and drag BellaBambi from the side of your nose to your ear. //

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Place BellaBambi in the middle of your chin under your bottom lip and drag it over your jaw to your earlobe. Repeat until you've finished your chin. Finally place BellaBambi below your jawbone and drag it from your upper neck to your collarbone. Please do not treat your larynx.

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Place BellaBambi below your toes and drag it to your ankle. You can also set it between your toes' extensions. Repeat on the bottom of your foot and drag BellaBambi to your heel.

Further applications on
www.bellabambi.com